



SurgiTrack

User manual

Surgical instrument and Set traceability

Version 0.29.6 · 09/10/2026 · www.surgitrack.eu

Contents

1. Getting started
2. From purchase to go-live
3. Set / instrument card
4. Overview
5. Department Sets & Instruments
6. Sterilization
7. Sterile expiry
8. Recycle bin
9. Connected devices
10. Instruments
11. Instrument Sets
12. Bulk import
13. Name check
14. Standalone in Use
15. Instrument Stock
16. Replacements & orders
17. Issues
18. Reports & Printing
19. Traceability
20. History
21. Management Center
22. Users & departments
- Glossary

Roles

Hospital admin	Users, departments, Studio and every screen of the hospital.
Sterilization supervisor	Everything a Sterilization user does, plus registering Sets and instruments, composition, Service.
Sterilization user	Receipt, check and composition, sterilizer, release, delivery.
Department user	Their department's Sets and instruments: dispatch, receipt, count, problem reports.
Viewer	Sees the whole hospital without changing anything.

Getting started

Your first sign-in, what you see on screen and where to ask for help.

Who it is for: All users

Step by step

1. Open the approval email and keep your username.
2. Press "Set password" and choose a password.
3. Sign in and open Help (the book icon, top right) on any screen you do not know.

Account

Only your hospital administrator opens accounts. If you got an invitation or a signup link, you fill in your name and department and wait for approval; you do not set a password yet.

Several devices

You can work from several computers and tablets at once. Each screen takes what the others saved or deleted at once (live sync); when the live connection is not available it checks every 20 seconds and as soon as you come back to the window. Your own changes are saved first. If two devices change the same Set or instrument at the same time, both changes are kept when they concern different fields; when they concern the same field, the one saved first is kept and a "Changed on another device at the same time" message names the field. If the network drops, a "No network connection" notice appears at the bottom of the screen: what you record is kept on this device and saved as soon as the connection is back, even if the page closes (open it again with the same user). The device keeps a copy of the hospital's data, so the app opens at once and fetches only what changed; the copy is wiped at sign-out and not used when older than a day. The top bar shows the state: "Live" (green) or "Sync" (grey, checking every 20 seconds), "Saving · N changes" while saving and "Not saved · retrying" if something failed; hovering over it shows the time of the last contact with the server. On a tablet the lists keep their key columns: long names take two lines, an issue's note shows under the item and, in portrait, each History entry shows the action and the user on a second line.

First sign-in

Once approved you get an "Access approved" email with your username (e.g. GN1234), your role and your department. Press "Set password", choose a password of at least 8 characters and sign in with your username or email.

Your screen

The menu on the left shows only what your role allows. After signing in the Overview opens if your menu has it; otherwise the first entry of your menu. In the scan field of the top bar, scan or type a barcode and the Set or instrument opens at once. Top right: language (EL/EN), Help, accessibility (text size, contrast) and notifications.

Unsaved changes

If you leave a page with details you have not saved (a new Set or instrument, editing a card's details, role permissions), the app first asks "Some changes are not saved". "Stay" is the default; "Leave without saving" discards them. The same goes for signing out; if you close or reload the tab, the browser asks.

Forgotten password

On the sign-in page press "Forgot password" and enter your email; you will get a link to a new password. After many wrong attempts sign-in locks for 15 minutes. The administrator can also send you a password change link.

Checks

- Do you see the right department and role? If not, tell the hospital administrator.
- Do not share your password: every movement is recorded under your name.

Tip. Help always opens on the page of the screen you are on. For access, password or permissions ask the hospital administrator; for a technical problem see "Support".

From purchase to go-live

The hospital administrator's steps, in order, from the first account to the first day of real use.

Who it is for: Hospital administrator

Step by step

1. Sign in as administrator from the Exeltos invitation.
2. Departments, libraries and sterilization flow.
3. Bulk import of instruments and Sets, then Name check.
4. Staff invitations and approvals.
5. Printer, scanners, labels.
6. Training in the Demo and a start day.

Your account

Exeltos creates your hospital and sends you an administrator invitation. From the email you set a password and sign in. If the hospital is on a trial, the app locks when it ends until activation.

Departments and libraries

In Users & departments → Departments add every department (Operating theatre, ICU, Central Sterilization...). In Studio → Libraries check specialties, manufacturers, sterilizers and color tapes; in Sterilization Flow the stages.

Instruments and Sets

Prepare the Excel file with the "Excel template" of Bulk import (Instruments or Sets → "Bulk import"): one row per instrument, the Set name on each of its rows. Check the preview and import. Then run the "Name check" so each code has one name.

Staff

In Users & departments → Users → "Add user": "Invite a person" one by one (email or link), "Hospital signup link" for many, or "From a CSV file". First the Sterilization supervisor, then Sterilization and the departments. Approve the requests (orange) with the right role and department.

Equipment

A barcode label printer in Sterilization, a barcode scanner (USB or Bluetooth, works as a keyboard) at every handover point, a computer or tablet with Chrome or Edge. If sterilizers will be connected, see "Connected devices". For the IT department: the app talks only to its own site and to Supabase (https and wss on *.supabase.co), loads nothing from third parties and does not open inside other sites; the site's security headers enforce this.

Training and go-live

Train staff in the Demo: the same screens with test data; you turn Demo access on per user in Users & departments. Label the Sets and instruments, set a start day and from that day every handover and receipt is done by scanning.

Checks

- Every department exists and is active.
- Sets have the right composition and department; the Name check shows nothing pending.
- Every user has a role and department; at least one Sterilization supervisor.
- The printer prints a label and the scanner opens the Set from the top bar.
- One test dispatch and receipt was completed by a department and Sterilization.

Tip. For help with any step, write to Exeltos at info@exeltos.com (see "Support").

Set / instrument card

Everything about one Set or instrument: details, composition, history, photos, uses and actions.

Who it is for: All users (actions depend on the role)

Step by step

1. Check the state and remaining uses.
2. See the history for the last movement.
3. Use the action you need from the top right.

SET CARD
ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ
ORTHO-BASIC

At department

DETAILS [Edit](#)

Status: At department

Sterilized / expires: 06/08/2026 06/02/2027

Code: ORTHO-BASIC

Name: ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ

Specialty: Orthopaedics

Department: Operating Theatre

Manufacturer: AESCULAP

Ownership: —

Color marker: No color

Usage type: Limited-use

Initial usage limit: 60

Uses: 24

Uses left: 36

Set barcode: 5000321

Set contents: 14/15 instruments

No photo
Capture or upload

Composition **History** **Issues 2** **Summary** **Notes**

SET COMPOSITION
Physical instruments
Each barcode is a separate physical instrument. The list uses all the available space.

+ Add instruments

14/15 Missing 1

Instrument, code, barcode or serial...

Filters

Instrument	Code	Name	Uses	Report
1	T001201	KOCHER ΕΥΘΕΙΑ 20 CM DEWIMED · S/N SN-20001	10/50	
2	T001202	ΨΑΛΙΔΙ ΜΑΥΟ ΕΥΘΥ 18 CM AESCULAP · S/N SN-20002	11/50	Open report
3	T001203	ΒΕΛΟΝΟΚΑΤΟΧΟ ΜΑΥΟ-HEGAR 20 CM DEWIMED · S/N SN-20003	12/50	
4	T001204	ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM AESCULAP · S/N SN-20004	13 uses	
5	T001205	ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM KLS MARTIN · S/N SN-20005	14 uses	Open report
6	T001206	ΑΓΚΙΣΤΡΟ LANGENBECK AESCULAP · S/N SN-20006	15 uses	
7	T001207	ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΛΑΒΙΔΑ 5MM KARL STORZ · S/N SN-20007	16/30	
8	T001208	KOCHER ΕΥΘΕΙΑ 20 CM DEWIMED · S/N SN-20008	17/50	
9	T001209	ΨΑΛΙΔΙ ΜΑΥΟ ΕΥΘΥ 18 CM AESCULAP · S/N SN-20009	18/50	
10	T001210	ΒΕΛΟΝΟΚΑΤΟΧΟ ΜΑΥΟ-HEGAR 20 CM DEWIMED · S/N SN-20010	19/50	
11	T001211	ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM AESCULAP · S/N SN-20011	20 uses	

Details

On the left: barcode, code, department, specialty, manufacturer, ownership and color marker. Users with the right can change them with "Edit".

Uses and limits

Limited-use instruments show their remaining uses. Each dispatch after use takes one use; at zero the instrument is taken out of use automatically.

History and photos

The tabs show every movement, problem report and photo of the item.

Actions

Depending on the role: send to Sterilization, print the barcode and, under "Actions", report a problem and manage (department, Set, Stock, Service, loss). A management action can be undone from its notice for a few seconds; the history keeps the undo too.

Previous and next card

When you open a card from a list (Sets, Instruments, Standalone, Stock), "3 of 24" with < > arrows shows next to "Back to list". You move to the previous or next card in the list's order and filters, also with Alt+← / Alt+→ on the keyboard. "Back to list" always returns to the list.


Overview

The whole hospital at a glance: Sets, instruments, sterilization flow, issues and usage limits.

Who it is for: Hospital administrator and Sterilization supervisor

Step by step

1. Start with the orange indicators: they need attention.
2. Click an indicator to see exactly which records it counts.
3. Check the latest movements for recent activity.

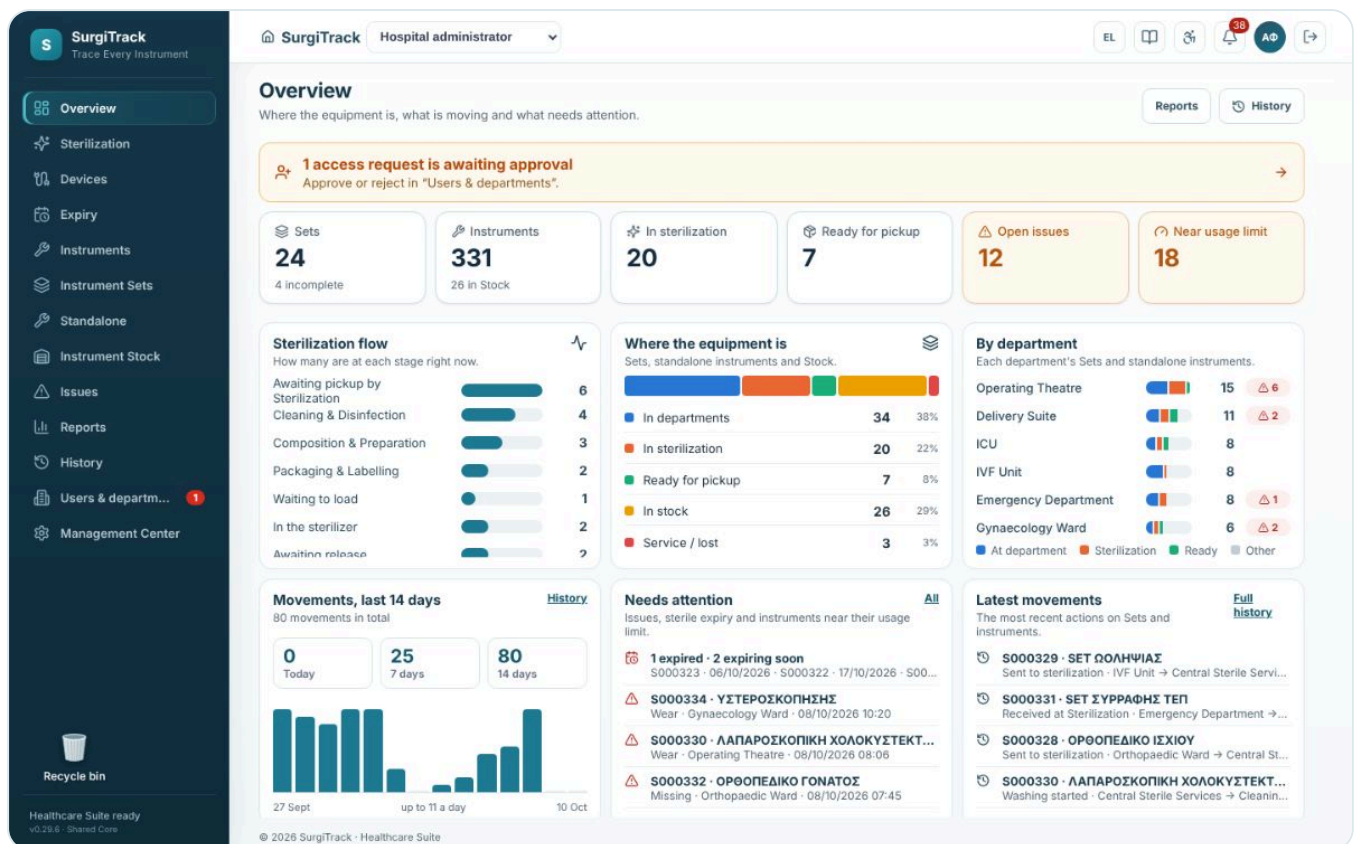

Indicators (KPIs)

Every indicator at the top is a link. It opens the matching list with its filters already applied, e.g. "Near usage limit" opens Instruments filtered to "Few uses left".

By department

The table shows what each department holds: at the department, in sterilization, ready and issues. Click a department name to see its Sets.

Sterilization flow

The bars show how many items are at each stage. Clicking a stage opens Sterilization on that tab.

Tip. The Overview is where the day starts: anything red or orange needs action today.

Department Sets & Instruments

The department's Sets and standalone instruments, and sending them electronically to Central Sterilization.

Who it is for: Department users

Step by step

1. Find the Set or instrument (search or scan the barcode).
2. Press "Send to Sterilization".
3. Enter the patient code where asked and confirm.
4. When it shows "Ready for pickup", collect it from Sterilization, signing with your user code and password.

The screenshot displays the SurgiTrack 'Operating Theatre' interface. The top navigation bar includes the SurgiTrack logo, a home icon, and user information (Nikos Δημητρίου). The main header shows the department's sets and instruments, and sending them to Sterilization. Below this, there are four status filters: 'At department' (8), 'In sterilization' (6), 'Ready for pickup' (1), and 'Open issues' (6). The interface is divided into two tabs: 'Instrument set' (8 Department sets) and 'Standalone instruments' (7 instruments in standalone use). A search bar is present above the list of instrument sets. The list includes details for various sets, such as 'ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ' (ORTHO-BASIC - Orthopaedics - 14/15 instruments), 'ΛΑΠΑΡΟΣΚΟΠΙΚΟ ΒΑΣΙΚΟ' (LAP-GEN - General Surgery - 12/12 instruments), 'ΝΕΥΡΟΧΕΙΡΟΥΡΓΙΚΟ ΒΑΣΙΚΟ' (NEURO-01 - Neurosurgery - 10/11 instruments), 'ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΧΟΛΟΚΥΣΤΕΚΤΟΜΗ' (LAP-CHOLE - General Surgery - 13/13 instruments), 'ΟΥΡΟΛΟΓΙΚΟ TUR' (URO-TUR - Urology - 11/11 instruments), and 'ΓΕΝΙΚΗΣ ΛΑΠΑΡΟΤΟΜΙΑΣ' (GEN-LAPAROT - General Surgery - 18/18 instruments). Each set entry shows its status (e.g., 'At department', 'Ready for pickup', 'In preparation', 'Cleaning & Disinfection', 'Packaging & Labelling', 'In sterilization') and a count of instruments used (e.g., 24/60, 17/50, 0 uses).

What I see

The two tabs list your Sets and standalone instruments. The indicators above filter the list: "At department", "In sterilization" (on the way to or inside Sterilization), "Ready for pickup".

Surgical count

In operating theatres (Studio → Libraries → Departments → "Count") "Send to Sterilization" asks for the count first. The "Count" button opens the form: Sterilization's part on top (composition, release, sterilization and expiry, cycle), the instruments below with "Sent" and "Counted". Enter the patient code, tick or scan each instrument or press "All present", then "Sign the count" (name and time).

Sending is then enabled. If instruments are missing, an issue opens and Sterilization sees it at Receipt. The form prints blank or signed, also from the Set's Print menu.

Send to Sterilization

After use, press "Send to Sterilization". The patient code is optional unless the Set or instrument has limited uses.

Instruments with a usage limit

For limited multi-use instruments (e.g. robotic), sending requires a patient code and confirming it was used: one use is recorded and the uses left are shown. When no uses are left the instrument is taken out of use automatically.

Report a problem

If an instrument is missing or damaged, open the Set, press "Actions" → "Report a problem" and describe it. Sterilization sees it under Issues.

Checks

- The patient code contains no patient name.
- The use was recorded for limited-use instruments.

Tip. Scan the barcode instead of searching by name: it is faster and avoids mistakes.

Sterilization

The Central Sterilization flow: receipt, washing, preparation, packaging, sterilizer load, load release, storage and delivery.

Who it is for: Sterilization users and supervisor

Step by step

1. Open the tab of your stage.
2. Scan or select the item.
3. Complete the stage check and move it on to the next stage.

The screenshot shows the SurgiTrack interface for a hospital administrator. The left sidebar contains navigation options: Overview, Sterilization (selected), Devices, Expiry, Instruments, Instrument Sets, Standalone, Instrument Stock, Issues, Reports, History, Users & departm..., and Management Center. The main area is titled 'Sterilization' and shows 'Your hospital's flow from Studio, fully traceable.' Below this is a 'Search & scan' bar with a barcode icon and a search input field. To the right of the search bar are buttons for 'Open' and 'Filters'. Below the search bar is a horizontal bar with stage tabs: Receipt (6), Cleaning (4), Composition (3), Packaging (2), Load (1), In the sterilizer (2), Release (2), and Delivery (7). Each tab has a sub-label: 'Chain of custody', 'Washing / disinfection', 'Instruments / deviations', 'Barrier / labelling', 'Sterilizer / load', 'Cycle running', 'Quality gate', and 'To department'. Below the tabs is a section titled 'Awaiting physical receipt' with 6 records. A 'Batch receipt' button is in the top right of this section. The table below lists the items:

Item	Department	Composition / status	Actions
S000327 SET SET ΕΠΙΣΙΟΤΟΜΗΣ	Delivery Suite Obstetrics	Composition 8 / 8	Receipt
S000328 SET ΟΡΘΟΠΕΔΙΚΟ ΙΣΧΙΟΥ	Orthopaedic Ward Orthopaedics	Composition 16 / 16 Count: 1 missing	Receipt
S000329 SET SET ΩΟΛΗΨΙΑΣ	IVF Unit Gynaecology	Composition 8 / 8	Receipt
S000344 SET ΣΕΤ ΑΝΑΣΘΕΩΦΟΝΗΣΗΣ ΠΑΙΔΙΚΟ (AMBU)	ICU Anaesthesiology	Composition 4 / 4	Receipt
T001311 INSTRUMENT ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM	Delivery Suite General Surgery	Standalone instrument	Receipt
T001486 INSTRUMENT KOCHER ΕΥΘΕΙΑ 20 CM	ICU General Surgery	Standalone instrument	Receipt

At the bottom left, there is a 'Recycle bin' icon and text: 'Healthcare Suite ready v0.29.6 Shared Core'. At the bottom center, it says '© 2026 SurgiTrack · Healthcare Suite'.

Stage tabs

Each tab is a stage and its number shows how many are waiting. Stages disabled in Studio are not shown. The blue "Search & scan" bar filters the tab as you type; scanning a barcode and pressing Enter opens the item in its stage. The Filters (department, specialty, Sets or instruments) sit next to it.

Physical receipt

In "Receipt" scan the Set and declare any visible deviation. You are recorded automatically as the receiver; the department person handing it over signs with their user code and password. Deviations (missing, damaged) create an issue.

Issues during Composition

Whatever needs action shows as a coloured row in the Set's instrument list, with a "Handle" button. An instrument with damage or another issue: replace it, send it to Service or Stock, or move it to another Set. A missing instrument ("Missing" row): "Cover the shortage" with an instrument from Stock, another Set or a standalone one, "Report", or "Keep as is" so the Set goes on with fewer instruments (the same as "Accept the recorded shortage" in the composition notice). An issue on the Set itself: "Fixed" with a note for the history, or "Keep as is" (the issue stays open).

Delivery to the department

In "Delivery" you are recorded automatically as the person handing over. The department person collecting signs with their user code and password on the same screen; without the signature the delivery cannot be completed.

Sterile shelf life

At "Packaging & Labelling" you choose a shelf life of 2, 3 or 6 months for each Set or instrument; the hospital default is suggested. It counts from the release and the expiry date shows on its card. If the hospital has no packaging stage, the shelf life is chosen at preparation.

Sterilizer load

In the "Load" tab every ready Set and instrument is already ticked; untick what does not go in and press one button, "Sterilizer load". Pick the sterilizer from the hospital list, the cycle number and the program, and declare which indicators are in the load (chemical, biological). No indicator result is given here. "Start cycle" puts the load in the sterilizer. If you pick a cycle from a connected device, that cycle has already ended and the load goes straight to Release; if the device reported a failure, it goes back to reprocessing.

In the sterilizer

The "In the sterilizer" tab shows the loads being sterilized now: sterilizer, cycle, program, who loaded it and when, and the items. They are locked; nobody can change them. When the cycle ends press "End of cycle": the load moves to Release and the screen goes there (Release is always required); "Cycle failed" sends the whole load back to reprocessing.

Load release

In the "Release" tab one button, "Release load", opens the loads that are waiting; pick the load and sterilizer. Confirm the physical parameters and the packaging and record the result of the indicators placed in the load. One passed indicator is enough, unless the unit requires a specific one. A failed indicator or "Do not release" sends the whole load back to reprocessing. Closing the release moves the items to "Delivery" to be handed to the department (or first to Storage, if the hospital uses it).

Release form

The "Form" button (in "In the sterilizer", in Release and in the recently released loads) prints the load's A4 form: load type, sterilizer, cycle, program, shelf life, the indicators with a place to stick their strip, the checks, the items with their expiry, the decision and the approval with the user's details and signature. From "In the sterilizer" it prints blank to fill in; from Release it prints filled in with what you have ticked there (indicators, checks) and the signed-in user's details; after the release it prints with

the recorded decision. Items show the sterilization date (STERILE) and the expiry date (hourglass); when no shelf life is set, the hospital default is shown.

Printing

The composition sheet (A4, with a check column and signatures) and the barcode label print from the preparation step. The label comes in three sizes (50×25, 70×35, 100×50 mm) with a SurgiTrack, logo, custom text or no header; the administrator sets the hospital default.

Checks

- The Set composition was checked before packaging.
- The load was recorded with sterilizer, cycle number and indicators.
- At release at least one indicator was recorded as passed.

Tip. An item under an active recall or with no uses left cannot circulate: the app stops it.

Sterile expiry

Sterile Sets and instruments with their expiry date: the ones expiring soon and the expired ones.

Who it is for: Sterilization and Administrator

Step by step

1. Press "Expiring soon" or "Expired" to see only those.
2. Filter by department and open the Set or instrument.
3. Ask the department to send it to Sterilization before it is used.

SurgiTrack Hospital administrator

Steriile expiry
Sets and instruments expiring or expired; once expired they need reprocessing.

Need attention: 3 | Expiring soon: 2 | Expired: 1 | Sterile: 5

Search: Name or barcode... Filters

Name	Barcode	Department	Location	Length	STERILE	Sterilization	Expires	Remaining
ΓΥΝΑΙΚΟΛΟΓΙΚΗΣ ΛΑΠΑΡΟΤΟΜΗΣ Set	S000323	Gynaecology Ward	At department	3 months	06/07/2026	06/10/2026	Expired 4 d ago	>
ΛΑΠΑΡΟΣΚΟΠΙΚΟ ΒΑΣΙΚΟ Set	S000322	Operating Theatre	Ready for pickup	2 months	17/08/2026	17/10/2026	Expires in 7 d	>
SET ΦΥΣΙΟΛΟΓΙΚΟΥ ΤΟΚΕΤΟΥ Set	S000325	Delivery Suite	At department	6 months	03/05/2026	02/11/2026	Expires in 23 d	>

Recycle bin

Healthcare Suite ready v0.29.6 - Shared Code

© 2026 SurgiTrack - Healthcare Suite

When it warns

A Set or instrument is "expiring soon" in the last month of its sterility (the last 10 days for 2 months). It then shows in the bell, in the Overview "Needs attention" and in this list; the department sees it in its own notifications.

When it expires

The app only warns and keeps it in the "Expired" list; it does not send it for reprocessing by itself. Once it is sent to Sterilization again, the expiry is cleared and a new one is set at the next release.

Tip. The hospital's default shelf life is set in Studio → Settings.


Recycle bin

Whatever is deleted stays 30 days in the bin and can be put back as it was.

Who it is for: Sterilization and Administrator

Step by step

1. Press the bin at the bottom left.
2. Find the record and press "Restore".

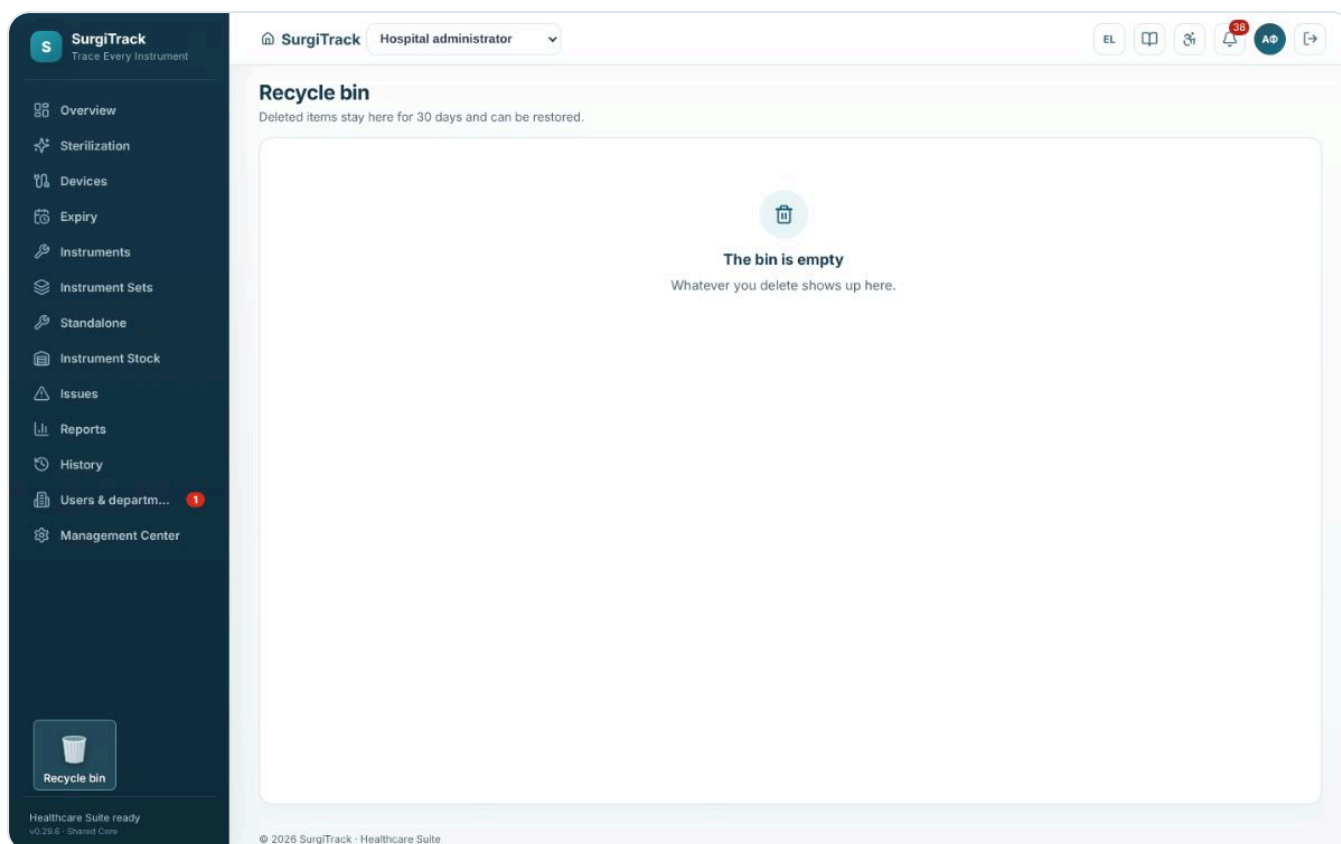

Where it is

The bin sits at the bottom left of the menu, like on a computer. When something is in it, it shows how many.

Restore

Sets, instruments, library records, color tapes and devices come back as they were; after 30 days they are deleted for good. Only the Administrator restores library records.

Connected devices

Sterilizers, washers and other devices that send their cycle data to SurgiTrack.

Who it is for: Sterilization and Administrator

Step by step

1. The admin adds the device and picks how it connects.
2. For network: "Network key", given to the manufacturer's technician.
3. For file or cable: "File" or "Cable" on the device card.

Connected devices
Sterilizers, washers and other devices that send cycle data to Sterilization.

STERILIZER Κλίβανος 1
MATACHANA · S1000 · S/N MT-2291
Network (automatic) · ...k3f9
Last cycle 1042 · 09/10/2026 21:39

WASHER-DISINFECTOR Πλυντήριο 2
MIELE · PG 8528 · S/N MI-7781
File (USB)
No data received yet

Cycles · Κλίβανος 1
"Recorded" means the cycle is already used in a Sterilization record.

Start	Cycle	Program	Result	°C	bar	Min	Source	Use
09/10/2026 21:39	1042	134°C 5'	Pass	134.6	3.1	52	Network (automatic)	Unused

Three ways to connect

Network: the device or the manufacturer software sends every cycle automatically with a device key.
File: export the cycles to USB (CSV or Excel) and upload them. **Cable:** the device connects to a Sterilization computer (Chrome or Edge) and what it prints is read automatically.

Cycle auto-fill

When recording a cycle or a load, the devices' new cycles are listed. "Use" fills in sterilizer, cycle number and program; you check and confirm. A cycle the device failed is marked as failed.

Record

Each device keeps its cycles with temperature, pressure, duration and result. "Download" saves them as Excel or CSV for review or audit.

Instruments

The register of every physical instrument, in Sets, standalone or in Stock.

Who it is for: Sterilization and Administrator

Step by step

1. Search by name, code, barcode or serial.
2. Narrow down with the filters.
3. Open the instrument for details and actions.

SurgiTrack Hospital administrator

Instruments
Every physical instrument: in stock, in a set or on its own.

Total instruments: **331**
 In sets: **265**
 Standalone in use: **40**
 Stock: **26**

Q Name, code, barcode, serial... Filters

Name	Barcode	Specialty	Location	Uses	Status
KOCHER ΕΥΘΕΙΑ 20 CM DEWIMED · 12.320.20 · S/N SN-20001	T001201	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	40 / 50	At department
ΨΑΛΙΔΙ MAYO ΕΥΘΥ 18 CM AESCULAP · 08.280.18 · S/N SN-20002	T001202	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	39 / 50	At department
ΒΕΛΟΝΟΚΑΤΟΧΟ MAYO-HEGAR 20 CM DEWIMED · 24.180.20 · S/N SN-20003	T001203	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	38 / 50	At department
ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM AESCULAP · 10.110.16 · S/N SN-20004	T001204	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	No limit	At department
ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM KLS MARTIN · 11.410.16 · S/N SN-20005	T001205	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	No limit	At department
ΑΓΚΙΣΤΡΟ LANGENBECK AESCULAP · 42.300.01 · S/N SN-20006	T001206	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	No limit	At department
ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΛΑΒΙΔΑ 5MM KARL STORZ · 70.510.05 · S/N SN-20007	T001207	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	14 / 30	At department
KOCHER ΕΥΘΕΙΑ 20 CM DEWIMED · 12.320.20 · S/N SN-20008	T001208	General Surgery	Set S000321 · ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ	33 / 50	At department

Recycle bin

Healthcare Suite ready
v0.29.6 · Shared Code

© 2026 SurgiTrack · Healthcare Suite

Filters

All filters are in the "Filters" button: department, specialty, manufacturer, location, state and usage type. The indicators above the list filter with one click.

Usage type

"With a usage limit" shows instruments with a usage limit, "Few uses left" those near it and "No limit" the rest.

Instrument card

Opening an instrument shows history, photos, uses, color marker and ownership. The limit and uses can be changed only by the Administrator and the Sterilization supervisor.

Import and names

"Bulk import" brings instruments and Sets from Excel. "Name check" finds names written differently and codes with many names, and unifies them.

Tip. The list loads progressively as you scroll; to find something quickly, use search.

Instrument Sets

The register of Sets with composition, department, state and completeness.

Who it is for: Sterilization and Administrator

Step by step

1. Find the Set by search or filters.
2. Open it for composition, history and actions.
3. Print the composition sheet when needed.

Instrument set
Set registry with separate Name, Code and unique Barcode fields.

Total Sets: 24 | At department: 7 | Ready for pickup: 4 | Sets with shortage: 4

Set name, code or barcode...

Set name	Barcode	Specialty	Department	Instrum...	Status	Expiry
ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ AESCULAP - ORTHO-BASIC	S000321	Orthopaedics	Operating Theatre	14 / 15	At department	06/02/2027
ΛΑΠΑΡΟΣΚΟΠΙΚΟ ΒΑΣΙΚΟ KARL STORZ - LAP-GEN	S000322	General Surgery	Operating Theatre	12 / 12	Ready for pickup	17/10/2026
ΓΥΝΑΙΚΟΛΟΓΙΚΗΣ ΛΑΠΑΡΟΤΟΜΗΣ DEWIMED - GYN-LAP	S000323	Gynaecology	Gynaecology Ward	13 / 14	At department	06/10/2026
ΝΕΥΡΟΧΕΙΡΟΥΡΓΙΚΟ ΒΑΣΙΚΟ KLS MARTIN - NEURO-01	S000324	Neurosurgery	Operating Theatre	10 / 11	Composition & Preparation	—
SET ΦΥΣΙΟΛΟΓΙΚΟΥ ΤΟΚΕΤΟΥ AESCULAP - DELIVERY-01	S000325	Obstetrics	Delivery Suite	10 / 11	At department	02/11/2026
SET ΚΑΙΣΑΡΙΚΗΣ ΤΟΜΗΣ DEWIMED - CSECTION-01	S000326	Obstetrics	Delivery Suite	12 / 12	Ready for pickup	26/03/2027
SET ΕΠΙΣΙΟΤΟΜΗΣ KLS MARTIN - EPISIO-01	S000327	Obstetrics	Delivery Suite	8 / 8	Awaiting pickup by Sterilization	—
ΟΡΘΟΠΕΔΙΚΟ ΙΣΧΙΟΥ AESCULAP - ORTHO-HIP	S000328	Orthopaedics	Orthopaedic Ward	16 / 16	Awaiting pickup by Sterilization	—
SET ΩΛΑΗΦΙΑΣ KARL STORZ - IVF-OPU	S000329	Gynaecology	IVF Unit	8 / 8	Awaiting pickup by Sterilization	—
ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΧΟΛΟΚΥΣΤΕΚΤΟΜΗ KARL STORZ - LAP-CHOLE	S000330	General Surgery	Operating Theatre	13 / 13	Cleaning & Disinfection	—

Composition

Each Set has a template composition (what it should contain) and its actual instruments. When something is missing, the Set shows as "Missing items".

Managing composition

The Sterilization supervisor adds instruments from Stock or removes them. An instrument with its own color joining a Set asks whether it keeps its color.

Ownership

A Set can belong to the hospital, a doctor or someone else. The doctor's name is searchable.

Bulk import

Import instruments and Sets from Excel or CSV, with a preview before anything is written.

Who it is for: Administrator and Sterilization supervisor

Step by step

1. Instruments or Sets → "Bulk import".
2. Download the template and fill it in.
3. Drag the file onto the page.
4. Check the preview, problems and names.
5. Press "Import".

Bulk import of instruments and Sets
From Excel or CSV, one row per instrument. Checked before writing and undone as a whole.

1 File — **2 Columns** — **3 Check** — **4 Import**

Drop the file here or click to choose
Excel (.xlsx) or CSV - for "Γενικό Νοσοκομείο Δοκιμών"

No file ready?
Download the template and fill it in. Or upload your own: in the next step you say which column is what. [Excel template](#)

How the file is read
One row per instrument. For example:

SET	INSTRUMENT	DEPARTMENT	QTY
Laparotomy 1	Kocher forceps	Theatre	4
Laparotomy 1	Metzenbaum scissors	Theatre	2
—	Farabeuf retractor	ICU	1
—	Pean forceps	—	3

- Same Set name = one Set: here one Set with 6 instruments (4 + 2).
- No Set: a standalone instrument of its department, or Stock without one.
- Quantity 4 = four instruments, each with its own barcode. Missing barcodes are assigned for you.
- Two identical Sets: give each its own name or Set barcode.

Imports into this hospital
No imports yet.

The file

Download the "Excel template". One row per instrument: Set, code, instrument name (required), quantity, department, specialty, manufacturer, max uses, serial number. Same Set name = same Set; no Set = standalone instrument; no department = Stock.

Preview

Before importing you see how many Sets, standalone and Stock instruments will be created, the instrument list with search and the rows with problems. Columns are recognized automatically; if one is not found, you map it by hand.

Names

The import can write names uniformly (capitals without accents, units written the same way) and use the name the hospital already has for the same code, so no duplicate names are created.

Undo

Every import appears under "Imports in this hospital" and can be undone in full with "Undo import", as long as none of its items has been used.

Checks

- Department names match those in Users & departments.
- Rows with problems were fixed or are not needed.
- After importing: "Name check".

Name check

Consistent instrument names: spelling fixes and codes that carry more than one name.

Who it is for: Administrator and Sterilization supervisor

Step by step

1. Instruments → "Name check".
2. Apply the spelling fixes.
3. Go through the codes with many names; the instrument count opens their list.
4. Update the compositions if asked.

Name check
One name per instrument. Every change is logged and can be undone.

Instruments: **331** | Distinct names: **21** | Spelling fixes: **7** (50 instruments) | Codes with several names: **0** (0 with a likely wrong code)

Spelling fixes: **7** | Same code, several names: **0**

Safe fixes
Spacing, capitals without accents, letters typed on the wrong keyboard (e.g. a Latin H in a Greek word) and units (12cm → 12 CM). The meaning of the name does not change.

	NOW	BECOMES	INSTRUMENTS
☒	ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΛΑΒΙΔΑ 5MM	ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΛΑΒΙΔΑ 5 MM	42
☒	ΛΑΜΑ ΛΑΡΥΓΓΟΣΚΟΠΙΟΥ MACINTOSH No 3	ΛΑΜΑ ΛΑΡΥΓΓΟΣΚΟΠΙΟΥ MACINTOSH NO 3	2
☒	ΜΑΣΚΑ ΣΙΛΙΚΟΝΗΣ AMBU No 4	ΜΑΣΚΑ ΣΙΛΙΚΟΝΗΣ AMBU NO 4	2
☒	ΛΑΜΑ ΛΑΡΥΓΓΟΣΚΟΠΙΟΥ MACINTOSH No 4	ΛΑΜΑ ΛΑΡΥΓΓΟΣΚΟΠΙΟΥ MACINTOSH NO 4	1
☒	ΛΑΜΑ ΛΑΡΥΓΓΟΣΚΟΠΙΟΥ MILLER No 1	ΛΑΜΑ ΛΑΡΥΓΓΟΣΚΟΠΙΟΥ MILLER NO 1	1
☒	ΜΑΣΚΑ ΣΙΛΙΚΟΝΗΣ AMBU No 2	ΜΑΣΚΑ ΣΙΛΙΚΟΝΗΣ AMBU NO 2	1
☒	ΜΑΣΚΑ ΣΙΛΙΚΟΝΗΣ AMBU No 5	ΜΑΣΚΑ ΣΙΛΙΚΟΝΗΣ AMBU NO 5	1

Apply to 50 instruments

Recycle bin
Healthcare Suite ready v0.25.6 - Shared Core

Spelling fixes

Risk-free fixes: double spaces, accents, lower case, Latin letters inside a Greek word (and the reverse), "12cm" → "12 CM". Keep the ones you want ticked and press "Apply to ... instruments".

Same code

For each code with many names the app suggests the most used one. Choose which names join, the final name (or type another) and press "Unify". "Possibly wrong code" means another instrument with the wrong code: fix it from its card.

Set compositions

A rename also reaches the Set template compositions. If "Set compositions with old names" appears, press "Update compositions".

Tip. Every change can be undone from the notice that appears right after it.

Standalone in Use

Instruments used on their own in departments, outside any Set.

Who it is for: Sterilization and Administrator

Step by step

1. Filter by department or usage type.
2. Open the instrument for details.

The screenshot displays the SurgiTrack web application interface. On the left is a dark sidebar with navigation options: Overview, Sterilization, Devices, Expiry, Instruments, Instrument Sets, **Standalone** (highlighted), Instrument Stock, Issues, Reports, History, Users & departm..., and Management Center. The main content area is titled 'Standalone instruments in use' with a subtitle 'Instruments used on their own in departments, outside a set.' and a '+ New instrument' button. Below the title are four filter boxes: 'In use' (40), 'Departments' (8), 'With usage limit' (27), and 'Few uses left' (8). A search bar 'Q Name, code or barcode...' and a 'Filters' button are also present. The main table lists instruments with columns: Name, Barcode, Specialty, Department, Uses, and Status. The table contains 12 rows of instrument data, including details like 'KOCHER EYΘEIA 20 CM', 'ΨΑΛΙΔΙ MAYO EYΘY 18 CM', and 'ΒΕΛΟΝΟΚΑΤΟΧΟ MAYO-HEGAR 20 CM'. Each row shows the instrument name, barcode, specialty, department, usage count/limit, and status (e.g., 'At department').

Name	Barcode	Specialty	Department	Uses	Status
KOCHER EYΘEIA 20 CM DEWIMED - 12.320.20	T001294	General Surgery	ICU	38 / 50	At department
ΨΑΛΙΔΙ MAYO EYΘY 18 CM AESCULAP - 08.280.18	T001295	General Surgery	Emergency Department	37 / 50	At department
ΒΕΛΟΝΟΚΑΤΟΧΟ MAYO-HEGAR 20 CM DEWIMED - 24.180.20	T001296	General Surgery	Gynaecology Ward	36 / 50	At department
ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM AESCULAP - 10.110.16	T001297	General Surgery	Orthopaedic Ward	No limit	At department
ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM KLS MARTIN - 11.410.16	T001298	General Surgery	Operating Theatre	No limit	At department
ΑΓΚΙΣΤΡΟ LANGENBECK AESCULAP - 42.300.01	T001299	General Surgery	ICU	No limit	At department
ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΛΑΒΙΔΑ 5MM KARL STORZ - 70.510.05	T001300	General Surgery	Emergency Department	12 / 30	At department
KOCHER EYΘEIA 20 CM DEWIMED - 12.320.20	T001301	General Surgery	Gynaecology Ward	31 / 50	At department
ΨΑΛΙΔΙ MAYO EYΘY 18 CM AESCULAP - 08.280.18	T001302	General Surgery	Orthopaedic Ward	30 / 50	At department
ΒΕΛΟΝΟΚΑΤΟΧΟ MAYO-HEGAR 20 CM DEWIMED - 24.180.20	T001303	General Surgery	Operating Theatre	29 / 50	At department
ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΛΑΒΙΔΑ 5MM					

What it is

A standalone instrument belongs to a department and circulates on its own. Many robotic instruments with a usage limit are standalone.

Usage limits

The "With a usage limit" and "Few uses left" indicators filter the list. Remaining uses show on each row.

Instrument Stock

Available instruments outside Sets and departments, ready to complete a Set.

Who it is for: Sterilization

Step by step

1. Open "Sets with missing items" to see what is missing.
2. Open the Set and press "Add instruments".
3. Pick the instrument from the Stock tab.

Instrument stock
Physical instruments available for replacement, set assembly or standalone issue.

Stock instruments: 26 | Instrument types: 9 | With usage limit: 17 | Below minimum: 0

Stock instruments 20 | Minimum stock levels

Q Name, code, barcode or serial... Filters

Instrument	Barcode	Specialty	Manufacturer	Uses	Status
KOCHER ΕΥΘΕΙΑ 20 CM 12.320.20	T001280	General Surgery	DEWIMED	0 / 50	Stock
ΨΑΛΙΔΙ ΜΑΥΟ ΕΥΘΥ 18 CM 08.280.18	T001281	General Surgery	AESCLAP	0 / 50	Stock
ΒΕΛΟΝΟΚΑΤΟΧΟ ΜΑΥΟ-HEGAR 20 CM 24.180.20	T001282	General Surgery	DEWIMED	0 / 50	Stock
ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM 10.110.16	T001283	General Surgery	AESCLAP	—	Stock
ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM 11.410.16	T001284	General Surgery	KLS MARTIN	—	Stock
KOCHER ΕΥΘΕΙΑ 20 CM 12.320.20	T001285	General Surgery	DEWIMED	0 / 50	Stock
ΨΑΛΙΔΙ ΜΑΥΟ ΕΥΘΥ 18 CM 08.280.18	T001286	General Surgery	AESCLAP	0 / 50	Stock
ΒΕΛΟΝΟΚΑΤΟΧΟ ΜΑΥΟ-HEGAR 20 CM 24.180.20	T001287	General Surgery	DEWIMED	0 / 50	Stock
ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM 10.110.16	T001288	General Surgery	AESCLAP	—	Stock
ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM					

Recycle bin
Healthcare Suite ready v0.29.6 - Shared Code
© 2026 SurgiTrack - Healthcare Suite

Numbers and tabs

At the top are the Stock numbers (instruments, kinds, with a usage limit, below the minimum); tapping a card filters the list, and "Below the minimum" opens Minimum stock. Below them, the "Stock instruments" and "Minimum stock" tabs switch the list.

Add to a Set

Adding happens from the Set: open the Set, press "Add instruments" and pick from the Stock tab. The "Sets with missing items" indicator opens the Sets that need completing.

Replacements & orders

Instruments in Service, damaged or worn, lost or out of use: replace from Stock in one click when it holds the same instrument, otherwise record a purchase order.

Who it is for: Sterilization and Administrator

Step by step

1. Filter by reason, department or Stock availability.
2. Replace from Stock what it holds.
3. Record a purchase order for the rest and print it.
4. When they arrive, "Receive into Stock": the new instruments enter Stock with their own barcodes.

The screenshot shows the SurgiTrack interface for a Hospital administrator. The 'Issues' section is active, displaying a list of instruments that need replacement. The interface includes a sidebar with navigation options like Overview, Sterilization, Devices, Expiry, Instruments, Instrument Sets, Standalone, Instrument Stock, Issues, Reports, History, Users & departm..., and Management Center. The main content area shows a summary of issues: 7 Need replacement, 5 In Stock, 2 Not in Stock, and 1 On order. Below this is a table of instruments with columns for Instrument, Set / department, Reason, Stock, and Status. Each row includes a 'Replacement' button.

Instrument	Set / department	Reason	Stock	Status
ΑΙΜΟΣΤΑΤΙΚΗ ΛΑΒΙΔΑ PEAN 14 CM T001533 - 15.340.14 - KLS MARTIN	S000321 ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ Operating Theatre	Service	No	Needs replacement 01/10/2026 10:15
ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM T001534 - 10.110.16 - AESCULAP	S000323 ΓΥΝΑΙΚΟΛΟΓΙΚΗΣ ΛΑΠΑΡΟΤΟΜΗΣ Gynaecology Ward	Service	Yes - 5	Needs replacement 04/10/2026 12:40
ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM T001205 - 11.410.16 - KLS MARTIN	S000321 ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ Operating Theatre	Damage / wear Service	Yes - 3	Needs replacement 21/08/2026 11:15
ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM T001535 - 11.410.16 - KLS MARTIN	S000325 SET ΦΥΣΙΟΛΟΓΙΚΟΥ ΤΟΚΕΤΟΥ Delivery Suite	Loss	Yes - 3	Needs replacement 05/10/2026 09:05
ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΛΑΒΙΔΑ 5MM T001508 - 70.510.05 - KARL STORZ	— Delivery Suite	Out of use	Yes - 1	Needs replacement 13/09/2026 12:20
ΨΑΛΙΔΙ ΜΑΥΟ ΕΥΘΥ 18 CM T001202 - 08.280.18 - AESCULAP	S000321 ΟΡΘΟΠΕΔΙΚΟ ΒΑΣΙΚΟ Operating Theatre	Damage / wear Damage	Yes - 5	Needs replacement 06/10/2026 15:10
ΨΑΛΙΔΙ METZENBAUM ΚΥΡΤΟ 18 CM T001203 - 08.280.18 - AESCULAP	S000341 ΑΓΓΕΙΟΧΕΙΡΟΥΡΓΙΚΟ	Damage / wear	No	Needs replacement

What it lists

Instruments in Service, with an open damage or wear report, lost and out of use. For each: the Set it belongs to (or left), the reason, since when, and whether Stock holds the same instrument (same code or, with no code, the same name).

Replace from Stock

"Replace" puts a matching Stock instrument into the Set. If the damaged one was still in it, it leaves for Service. Select several and press "Replace from Stock" for all at once; one Stock instrument is never used twice.

Purchase order

For what Stock lacks: select it and "Purchase order". One line per kind with its quantity; add supplier and note. In the "Orders" tab the order prints, is marked ordered or cancelled, and "Receive into Stock" registers the new instruments in Stock.

Tip. "Download" (Excel, PDF / Print, CSV) takes your selection, otherwise what the filters show.

Issues

Missing, damaged, faulty or lost items that need handling.

Who it is for: All users

Step by step

1. Press the "Open" card (or "Open over 7 days" for those waiting long).
2. Open the issue and see its note and photos.
3. Resolve and close it.

SurgiTrack Hospital administrator

Issues
Reported problems and instruments that need replacing.

Problem reports 12 **Replacements & orders** 7

Total **16** Open **12** Open over 7 days **4** Resolved **4**

Instrument/set, barcode, type or note...

Item	Note	Created	Status
S000334 - ΥΣΤΕΡΟΣΚΟΠΗΣΗΣ Wear - Gynaecology Ward	Χαλαρή άρθρωση σε λαβίδα Kocher.	08/10/2026 10:20	Open
S000330 - ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΧΟΛΟΚΥΣΤΕΚΤΟΜΗ Wear - Operating Theatre	Φθαρμένη μόνωση σε λαπαροσκοπική λαβίδα 5mm.	08/10/2026 08:06	Open
S000332 - ΟΡΘΟΠΕΔΙΚΟ ΓΟΝΑΤΟΣ Missing - Orthopaedic Ward	Λείπει ένα άγκιστρο Langenbeck από τη σύνθεση.	08/10/2026 07:45	Open
T001531 - ΨΑΛΙΔΙ METZENBAUM ΚΥΡΤΟ 18 CM Wear - Operating Theatre	Φθαρμένες λεπίδες στην άκρη.	07/10/2026 11:30	Open
S000340 - SET ΜΙΚΡΟΕΠΕΜΒΑΣΕΩΝ ΤΕΠ Missing - Emergency Department	Λείπει ένα ψαλίδι ραμμάτων.	06/10/2026 18:45	Open
T001202 - ΨΑΛΙΔΙ MAYO ΕΥΘΥ 18 CM Damage - Operating Theatre	Χαλαρή άρθρωση, δεν κόβει καθαρά.	06/10/2026 15:10	Open
T001535 - ΛΑΒΙΔΑ ΧΕΙΡΟΥΡΓΙΚΗ 16 CM Loss - Delivery Suite	Δεν επέστρεψε με το Σετ μετά τον τοκετό.	05/10/2026 09:05	Open
T001534 - ΛΑΒΙΔΑ ΑΝΑΤΟΜΙΚΗ 16 CM Damage / Service - Gynaecology Ward	Αποστείρωση - σύνθεση & προετοιμασία: μεταφέρθηκε στα...	04/10/2026 12:40	Open
S000338 - SET ΤΡΑΧΕΙΟΣΤΟΜΙΑΣ Service - ICU	Επισκευή ψαλιδιού από τον προμηθευτή. Επέστρεψε.	04/10/2026 12:00	Resolved
S000341 - ΔΕΓΓΕΙΟΧΕΙΡΟΥΡΓΙΚΟ			

© 2026 SurgiTrack - Healthcare Suite

Creating

An issue is created from a Set or instrument card ("Actions" → "Report a problem") or automatically from a receipt deviation. Sterilization declares a loss or a Service send-off in "Actions" → "Manage", because they change the state; a department reports a loss and Sterilization confirms it.

Closing

Sterilization resolves the issue (replacement, service, return) and closes it. Everything is recorded in History.

Two tabs

Sterilization and the Administrator see two tabs: "Problem reports" (what was reported) and "Replacements & orders" (the instruments to replace, from Stock or with a purchase order).

Numbers

Above the list of reports, four cards: Total, Open, Open over 7 days (those waiting long) and Resolved. They count within the department, type and search you chose; pressing one shows only those in the list. The same date filter is also in "Filters".

Reports & Printing

Ready-made reports with filters and download as Excel, PDF or CSV.

Who it is for: Sterilization and Administrator

Step by step

1. Pick a report on the left.
2. Set the filters.
3. Press "Download" and choose Excel, PDF or CSV.

Available reports

Set composition, By department, By specialty, Service & faults, Usage limits, Out-of-use instruments and Patient traceability.

Download: Excel, PDF, CSV

The "Download" button has three options: "Excel" downloads the results as a table with filters, "PDF / Print" opens an A4 preview to print or save as PDF, and "CSV" gives them as plain text for other programs. All take the results with the filters you set. The table fits the screen width; long texts wrap.

Out of use

Instruments that reached their usage limit leave the lists and stay as history in this report, with date and reason.

Traceability

The full path of a Set or instrument, or everything used for one patient, on a timeline.

Who it is for: Sterilization, Supervisor, Admin, Viewer

Step by step

1. Scan the barcode or type the patient code.
2. Read the timeline from the most recent movement.
3. Download or print the results if needed.

SurgiTrack Trace Every Instrument

Hospital administrator

Tracing

Search by set/instrument barcode or patient code — without patient names.

Q S000330

Download

Search

SET
S000330
ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΧΟΛΟΚΥΣΤΕΚΤΟΜΗ · Operating Theatre · Cleaning & Disinfection

Open card

08/10/2026 08:05	● Washing started Central Sterile Services → Cleaning & Disinfection · Demo Sterilization User
08/10/2026 07:55	● Received at Sterilization Operating Theatre → Central Sterile Services · Demo Sterilization User
26/09/2026 13:57	● Deliver to department · receiver signature Central Sterile Services → Operating Theatre · Demo Sterilization User
25/09/2026 14:39	● Received at Sterilization Operating Theatre → Central Sterile Services · Demo Sterilization User
25/09/2026 13:03	● Sent to sterilization Operating Theatre → Central Sterile Services · Demo Operating Theatre User · Patient PT-2026-0129
18/09/2026 13:57	● Deliver to department · receiver signature Central Sterile Services → Operating Theatre · Demo Sterilization User
17/09/2026 09:39	● Received at Sterilization Operating Theatre → Central Sterile Services · Demo Sterilization User
17/09/2026 08:03	● Sent to sterilization Operating Theatre → Central Sterile Services · Demo Operating Theatre User · Patient PT-2026-0128
13/09/2026 13:57	● Deliver to department · receiver signature Central Sterile Services → Operating Theatre · Demo Sterilization User
12/09/2026 10:39	● Received at Sterilization

Recycle bin

Healthcare Suite ready
v0.29.6 · Shared Code

© 2026 SurgiTrack · Healthcare Suite

Search

Scan or type a Set (S...) or instrument (T...) barcode, or a patient code, and press "Search". Patients are searched by code only; no names are recorded.

Timeline

Each movement shows date and time, action, from where to where, user and, where there is one, the patient code. For a Set or instrument, "Open card" opens its card.

Download and print

The "Download" button gives the results as Excel, CSV or PDF for printing, e.g. for a quality review or a recall.

History

The chain of custody: every Set and instrument movement, who, when and from where to where.

Who it is for: All users (department users see only their department)

Step by step

1. Filter by date or barcode.
2. Click a movement for full details.

SurgiTrack Hospital administrator

Movement history
Unalterable history of handovers, receipts, checks and transfers.

Search: Barcode, set/instrument, user or patient code...

Showing the history of the last 90 days.

251 records

Date / time	Set / Instrument	Route	Action	User
08/10/2026 08:32	SET S000329 SET ΔΟΛΗΦΙΑΣ	IVF Unit > Central Sterile Services	Sent to sterilization	Demo IVF User Patient PT-2026-0191
08/10/2026 08:12	SET S000331 SET ΣΥΡΡΑΦΗΣ ΤΕΠ	Emergency Department > Central Sterile Services	Received at Sterilization	Demo Sterilization User
08/10/2026 08:10	SET S000328 ΟΡΘΟΠΕΔΙΚΟ ΓΥΝΑΙΚΕΙΟ	Orthopaedic Ward > Central Sterile Services	Sent to sterilization	Demo Χρήστης Ορθοπαιδικής Patient PT-2026-0188
08/10/2026 08:05	SET S000330 ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΧΟΛΟΚΥΣΤΕΚΤΟΜΗ	Central Sterile Services > Cleaning & Disinfection	Washing started	Demo Sterilization User
08/10/2026 07:55	SET S000330 ΛΑΠΑΡΟΣΚΟΠΙΚΗ ΧΟΛΟΚΥΣΤΕΚΤΟΜΗ	Operating Theatre > Central Sterile Services	Received at Sterilization	Demo Sterilization User
08/10/2026 07:40	SET S000332 ΟΡΘΟΠΕΔΙΚΟ ΓΟΝΑΤΟΣ	Cleaning & Disinfection > Inspection & Assembly	Washing completed	Demo Sterilization User
08/10/2026 07:20	SET S000333 ΟΥΡΟΛΟΓΙΚΟ TUR	Inspection & Assembly > Packaging & Labelling	Composition complete	Demo Sterilization User
08/10/2026 07:05	SET S000334 ΥΣΤΕΡΟΣΚΟΠΗΣΗΣ	Packaging & Labelling > Sterilization	Loaded into ΚΑΙΒΟΤΟ Αρμού 01	Demo Sterilization User
08/10/2026 07:05	SET S000335	Packaging & Labelling >	Loaded into ΚΑΙΒΟΤΟ Αρμού 01	Demo Sterilization User

© 2026 SurgiTrack - Healthcare Suite

Recent and older history

To open quickly, the app loads the last 90 days of history (compositions and releases for 7 months, as long as the longest sterile shelf life). In History, "Load older" brings the rest; Traceability and Reports load it by themselves. Nothing is deleted: all history stays in the database.

Who made the movement

Every entry keeps the account it was recorded under; the database itself writes it, not the device, and it never changes. The movement details show it as "Account (from the database)". When the recorded name is not the account's, the account also shows in the list, in orange.

Cleaning up history (owner)

Only the platform owner, inside a hospital, sees "Clean up history". They tick the entries (or "Select all shown" after filtering) and press "Delete selected". The deletion is permanent and is recorded as a new "History clean-up · N entries" entry. For everyone else the history never changes.

Date format

All dates show the same way, day/month/year and time in 24 hours: 08/10/2026 13:05. Entries recorded earlier as "8/10/26, 1:05 μ.μ." show this way too; the entry in the database does not change.

Immutable

History records are never changed. Use the date, route and action filters to find a movement.

Management Center

Hospital settings: libraries, workflow stages, role permissions, color tapes and the usage warning threshold.

Who it is for: Administrator

Step by step

1. Pick a section.
2. Make your changes.
3. Save; changes apply at once.

SurgiTrack Studio
Libraries, sterilization flow, role permissions and hospital settings.

LIBRARY Specialties
Surgical and clinical specialties.

Search name or code...

NAME	ENGLISH	CODE
Γενική Χειρουργική	General Surgery	-
Ορθοπαιδική	Orthopaedics	-
Γυναικολογική	Gynaecology	-
Νευροχειρουργική	Neurosurgery	-
Ουρολογία	Urology	-
Ανασθησιολογία	Anaesthesiology	-

Libraries

Departments, specialties, manufacturers, suppliers, instrument categories and sterilizers. Whatever you define here appears in the choices across the app. In "Sterilizers", "Add sterilizers" adds several at once in the hospital's naming (Sterilizer A, B, C... or 1, 2, 3...); the choice is kept for next time and the sterilizers appear in the Sterilizer load.

Permissions

Set what each role can do. Some permissions (creation, composition, usage limits) belong only to the Sterilization supervisor. The database itself also limits a department user: on Sets and instruments it changes only what its own workflow changes (dispatch to Sterilization from the department, uses, out of use at the limit, the count), and issues it can only report; anything else it sends stays as it was.

Sterilization Flow

The stages from receipt to release. Turn on only the stages and checks your hospital's Central Sterilization performs. In the release policy you set whether a specific indicator is required; without one, one passed indicator (chemical or biological) is enough.

Roles and Settings

"Roles" explains what each role does, so you choose correctly when approving. In "Settings" you set when an instrument shows as "Few uses left", the default sterile shelf life (2, 3 or 6 months) and the hospital's barcode label (size, header, logo), and "Screen lock after inactivity" (default 15 minutes, or Never): on shared tablets the screen locks when not in use and unlocks with the same user's password, or another user signs in as themselves; syncing goes on behind it and nothing is lost.

Users & departments

The hospital's users, invitations, access requests and departments.

Who it is for: Hospital administrator

Step by step

1. Departments: make sure they all exist.
2. Users → "Add user": invite a person, the hospital link or CSV.
3. Approve the requests with the right role.
4. Press "Refresh" for the latest requests.

Γενικό Νοσοκομείο Δοκιμών
The hospital's users, invitations and departments.

Users 9 1 Departments 10

Search name, email, department or role...

1 request awaits approval
Check the details, pick a role and department, then approve. The user already has a username and password; on approval they get one email and can sign in.

ΠΑΝΑΓΙΩΤΗΣ ΡΗΓΑΣ
prigas · p.prigas@hospital.test · asked for Operating Theatre · 09/10/2026 00:03

Role: Department user Department: Operating Theatre Note (optional):

Approve Reject

USER	USERNAME	DEPARTMENT	ROLE	ACCESS
Αριστέλης Φιλοκώστας admin@hospital.test	admin	Whole hospital	Hospital administrator	Active
Γεώργιος Αντωνίου supervisor@hospital.test	supervisor	Central Sterile Services	Sterilization supervisor	Active
Δημήτρης Λαμπρόπουλος staff3@hospital.test	staff3	Gynaecology Ward	Department user	Active
Ελένη Κωνσταντίνου viewer@hospital.test	viewer	Whole hospital	Viewer (read only)	Active
Κατερίνα Ιωάννου staff4@hospital.test	staff4	ICU	Department user	Inactive

© 2026 SurgiTrack · Healthcare Suite

Invitation

"Add user" → "Invite a person": role, email and (if you like) department. "Send email" sends the signup link; "Copy link" copies it to send any way you like (Viber, SMS). It is valid 7 days. The user fills in their name, department and password and sees their username at once. Until they sign up the invitation shows here, where you resend it, copy the link or cancel it.

Hospital link and CSV

For many at once: share the "Hospital signup link" (10 days) with staff; they sign up the same way and you approve them. "From a CSV file" (Name; Email; Department; Role) creates the accounts at once, without approval, and each person gets an email with a username and a password link.

Approval

Requests show in orange, with a count on the bell (no email is sent). Check the name, email and department, set the role and press "Approve" or "Reject". On approval the user gets one email and signs in with their username and the password they set; on rejection they get one email and their account is removed.

Users

Open a user to change the role, department or Demo access, or to deactivate them. "Send password change link" helps someone who forgot their password. "Delete" removes the account completely; their movement history stays.

Departments

Add, rename or deactivate departments. The names appear in signups, imports and every list; a deactivated department is not offered to new users.

Checks

- Every department user has the right department.
- People who left the hospital are inactive or deleted.

Tip. If an invitation email does not arrive, ask them to check spam or copy the link and send it yourself.

Glossary

Όνομα χρήστη / Username

Your personal sign-in code (e.g. GN1234: your initials and four digits). It comes with the approval email; you sign in with it or with your email.

Αίτηση πρόσβασης

A new user signup waiting for the hospital administrator to approve it. Until approval there is no account and no password.

Πρότυπη σύνθεση

What a Set should contain (code, name, quantity). It is compared with its actual instruments to show what is missing.

Έλεγχος ονομασιών

A tool that makes names consistent: the same spelling and one name per code.

Παρατηρητής / Viewer

A read-only account (e.g. Nursing Directorate): sees the overview, registries, issues, history and reports of the whole hospital and changes nothing. The hospital administrator creates it.

Σετ / Set

A group of instruments with a defined composition, sterilized and circulated together.

Μεμονωμένο / Standalone

An instrument that belongs to a department and circulates on its own.

Απόθεμα / Stock

Available instruments outside Sets and departments, used to complete Sets.

Όριο χρήσεων

How many times a limited multi-use instrument may be used. One use is counted on each dispatch after use.

Εκτός χρήσης

An instrument that reached its usage limit; it leaves the lists and stays as history in reports.

Κωδικός ασθενούς

An identifier for traceability, without the patient name.

Barcode

The unique scan code of each Set (S...) or instrument (T...).

Αποδέσμευση

The check after a sterilization cycle that allows the item to be used.

Φορτίο

The Sets and instruments that go into the same sterilizer cycle; they are released or sent back to reprocessing together.

Χημικός / βιολογικός δείκτης

Indicators placed in the load that show whether the cycle worked. They are declared at loading and their result is recorded at release; one passed indicator is enough.

Διάρκεια / λήξη αποστείρωσης

How long a Set or instrument stays sterile (2, 3 or 6 months from the release). It warns in the last month (10 days for 2 months).

Ανάκληση

A procedure that blocks items from a faulty cycle from circulating.

Χρωματική σήμανση

Colored tapes that identify a Set or instrument at a glance.